


# FUJINON VIDEO BRONCHOSCOPES *Designed to improve physician productivity and patient care.*

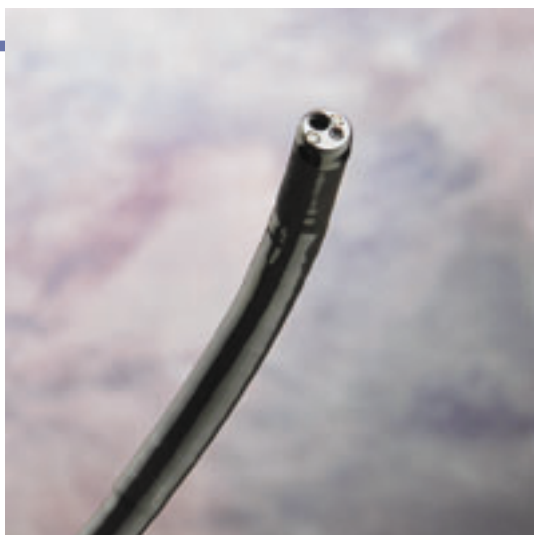

## EB-270S

The EB-270S sets a new standard in bronchoscopes with a slimmer distal end while maintaining a generous forceps channel. At the same time, high resolution video image provides clearer and more detailed images than ever before as a result of the advanced, miniature CCD chip built into the scope tip. Insulating polymer materials allow for full compatibility with high frequency treatment devices and YAG lasers.

## EB-270T

The EB-270T establishes a new benchmark in therapeutic bronchoscopes with a generous forceps channel within a slim outside diameter, allowing for greater therapeutic access. Suction performance has been enhanced and the forceps tube junction redesigned to facilitate even the most challenging procedures. Insulating polymer materials allow for full compatibility with high frequency treatment devices and YAG lasers.

	EB-270S	EB-270T
Viewing direction	Forward	
Observation range	3 ~ 100mm	
Field of view	120 °	
Distal end diameter	4.9mm	5.9mm
Flexible portion diameter	4.9mm	5.9mm
Bending capability : UP	180 °	
DOWN	130 °	
Forceps channel diameter	2.0mm	2.8mm
Working length	600mm	
Total length	870mm	

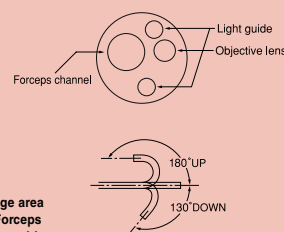


Image area & Forceps entry position



## EPX-2200 Xenon 150W Processor

### Specifications

TV system	NTSC / PAL
Power supply	AC120V/230V
Image zoom	Electronic Zoom X1.00 to X2.00 (0.05steps)
Memory	Patient data-44 patients Clinical procedure-20 procedures
Light source	150 watt Xenon lamp (Emergency lamp 75W Halogen lamp)
Dimensions	330(W) X 430(D) X 160(H)mm, 13 X 16.9 X 6.3 inches
Weight	16kg (35.2lbs)
Image output	Video(BNC) 1 channel S-Video 2 channels RGB-VIDEO 2 channels



### FUJI PHOTO OPTICAL CO., LTD.

1-324 UETAKE, KITA-KU, SAITAMA CITY, SAITAMA 331-9624, JAPAN  
TEL:81-48-668-2153 FAX:81-48-668-1570  
E-mail:sales@msv.fujinon.co.jp URL <http://www.fujinon.co.jp/>

### FUJINON INC.

10 HIGH POINT DRIVE, WAYNE, NJ 07470, U.S.A.  
TEL:1-973-633-5600 FAX:1-973-633-8818  
E-mail:med.info@fujinon.com URL <http://www.fujinon.com/>

### FUJINON (EUROPE) GmbH

HALSKESTRASSE 4, 47877 WILLICH, GERMANY  
TEL:49-2154-924-0 FAX:49-2154-924-290  
E-mail:Fujinon@fujinon.de URL <http://www.fujinon.de/>

### FUJI KOKI SINGAPORE PTE. LTD.

BLK 211 HENDERSON ROAD #10-04, HENDERSON INDUSTRIAL PARK,  
SINGAPORE 159552  
TEL:65-6276-4988 FAX:65-6276-6911 E-mail:fujinon@fujinon.com.sg

### F.P.FOCUSING PTY. LTD.

UNIT-18, 52 HOLKER STREET, SILVERWATER N.S.W.2128, AUSTRALIA  
TEL:61-2-9748-2744 FAX:61-2-9748-2428 E-mail:fpfocus@dot.net.au

### Shanghai FUJINON Medical Equipment Co. Ltd.

FLAT B, 6/FL, SHI MEI BLD, 79 RI JING ROAD,  
WAI GAO QIAO FREE TRADE ZONE, SHANGHAI 200131, CHINA  
TEL:86-21-5424-9262(Office) TEL:86-21-5866-2151(Service)  
FAX:86-21-5424-9264 E-mail:sh-fujinon@etang.com